


This fact sheet gives you important information about the general use and care of a swivel shower chair. It describes how to install, use and care for a swivel shower chair.

**This fact sheet does not replace specific instructions given to you by the equipment prescriber.**

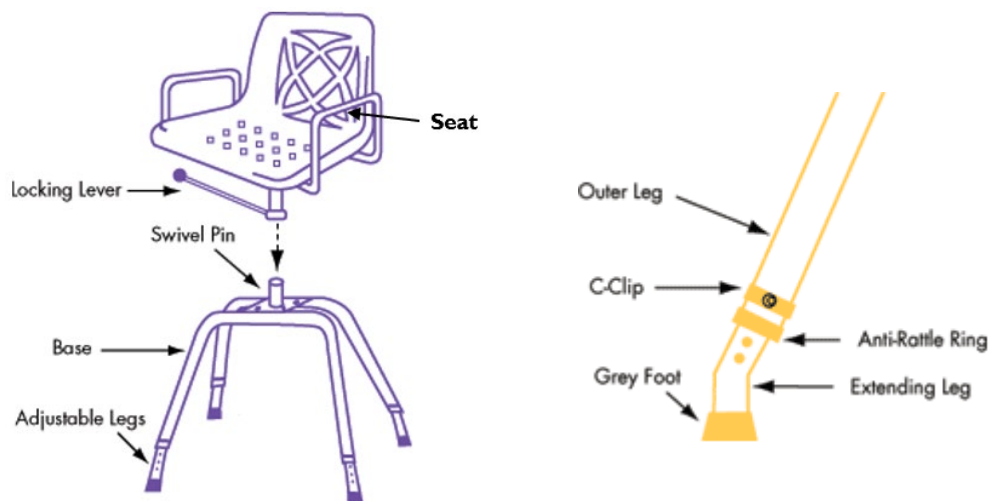

### Installing the swivel shower chair

1. Place the base of the swivel shower chair in the shower. Make sure the four non-slip feet are sitting flat on the ground.
2. Place the seat on the base.
3. There is a lever on the right hand side of the seat that locks the seat into place. To lock the seat into position, push the lever down to fit into the slot in the base.
4. To unlock the seat, lift the lever.

### Adjusting the height of the swivel shower chair

1. Adjust one leg at a time.
2. Remove the C-clip on the leg by pushing the open side and pulling out the pin.
3. Adjust to the height required – so that one of the holes in the extending leg is visible through the hole in the outer leg. Replace the C-clip. Ensure that the pin goes through both holes in the outer leg and the plastic ring is locked in position.

## Using the swivel shower chair

1. Turn the seat so that it faces the side of the shower and lock it into position.
2. Stand close to the shower, so you can feel the chair against the back of your legs.
3. Reach behind you and hold both of the armrests for support.
4. Lower yourself slowly into the chair using the armrests.
5. Unlock the seat by pulling up the lever on the right hand side.
6. Turn the seat and lift your legs into the shower.
7. Lock the seat in place by pushing the lever down.
8. Shower while sitting on the swivel chair.
9. To get out of the shower, reverse the process.

## Care of the swivel shower chair

- Wipe any soiling immediately with hot, soapy water.
- Dry the seat with a towel after each use.
- Clean the swivel shower chair regularly with hot, soapy water.
- Check for rust and signs of wear to the plastic seat.

## Repairs to your equipment

If your equipment is broken, damaged, rusted or worn contact your local TasEquip warehouse to arrange for repair or replacement. Please call 1300 827 378 or contact your regional warehouse directly:

### North

63 Dowling Street  
(entry off Henry Street)  
Launceston 7250  
Phone: (03) 6777 4336  
Email:  
[NthTasEquip@ths.tas.gov.au](mailto:NthTasEquip@ths.tas.gov.au)

### North West

Mersey Community Hospital  
Moriarty Road  
Latrobe 7307  
Phone: (03) 6478 5599  
Email:  
[NWTasEquip@ths.tas.gov.au](mailto:NWTasEquip@ths.tas.gov.au)

### South

Repat Centre  
90 Davey Street  
Hobart 7000  
Phone: (03) 6166 7393  
Email:  
[SouthTasEquip@ths.tas.gov.au](mailto:SouthTasEquip@ths.tas.gov.au)

Review date: September 2025